

UYA 4Y PLUS Ultra-Microbalances

MYA 4Y PLUS Microbalances



SYNERGY LAB Line – New quality of small mass weighing.



UYA 4Y PLUS



MYA 4Y PLUS



Automatically opened weighing chamber



Proximity sensors

Functions

Parts counting	Percent weighing	Statistical quality control	Proximity sensors	Multilingual menu
Dosing	Statistics	Autotest	Ambient conditions measurement	Replaceable unit
Checkweighing	Animal weighing	Automatic sliding door	GLP procedures	Moveable range
Formulations	Differential weighing			

Features

Excellent Readability Starting from 0.1 µg

Due to exceptional weighing parameters, the UYA 4Y PLUS and MYA 4Y PLUS microbalances are intended for the most demanding laboratory applications.

Intuitive Operation and Touch Screen

5.7" colour touch screen enables intuitive operation and easy access to numerous applications and functions of the weighing instrument.

Numerous Options of Data Management

Extensive storage capacity enables record of all measurement data in a form of complex reports and statistical graphs.

21 CFR Part 11 Compliance Module

MYA 4Y PLUS and UYA 4Y PLUS microbalances are compliant with 21 CFR part 11 and EU GMP part 4, annex 11. This means that the PLUS series features numerous functions for protecting and monitoring work results of particular operators.

Smart Min Weight

With use of the „Smart Min Weight“ function, the MYA 4Y PLUS microbalance and UYA 4Y PLUS ultra-microbalance automatically adjust resolution of the range to the currently weighed mass. This results with improvement of the minimum sample weight parameter by 30%.

REFLEX LEVEL SYSTEM

New autolevelling system ensures ground tilt test, automatic balance levelling, control of the correct balance level and confirmation of correct balance level on adjustment and weighing reports.

ALIBI Memory

Data security and protection is provided by ALIBI memory which automatically archives all carried out measurements.

Technical Specifications

	UYA 2.4Y PLUS	UYA 6.4Y PLUS	MYA 0,8/3.4Y PLUS	MYA 2.4Y PLUS
Maximum capacity [Max]	2.1 g	6,1 g	0.8 g / 3 g	2.1 g
Minimum load	10 µg	10 µg	100 µg	100 µg
Readability [d]	0.1 µg	0,1 µg	1 µg / 10 µg	1 µg
Verification scale interval [e]	1 mg	1 mg	1 mg	1 mg
Tare range	-2.1 g	-6,1 g	-3 g	-2.1 g
Standard values				
Repeatability (5% Max)*	0.15 µg	0,2 µg	0.6 µg	0.41 µg
Repeatability (Max)	0.35 µg	0,45 µg	4.1 µg	1 µg
Minimum weight (USP)	0.3 mg	0.4 mg	1.2 mg	0.82 mg
Minimum weight (U=1%, k=2)	0.03 mg	0.04 mg	0.12 mg	0.082 mg
Permissible values				
Repeatability (5% Max)*	0.35 µg	0,4 µg	1.2 µg	0.8 µg
Repeatability (Max)	0.6 µg	0,8 µg	6 µg	1.5 µg
Linearity	±1.5 µg	±1.5 µg	±3 µg / ±10 µg	±3 µg
Eccentric load deviation	1.5 µg	1.5 µg	3 µg / 10 µg	3 µg
Stabilization time	10 ÷ 20 s	10 ÷ 20 s	max 8 s	max 8 s
Adjustment	internal	internal	internal	internal
Moveable range	—	—	Yes	—
Verification	Yes	Yes	Yes	Yes
OIML Class	I	I	I	I
Indicator fastening	35 cm cable, wireless connection (option)**	35 cm cable, wireless connection (option)**	35 cm cable, wireless connection (option)**	35 cm cable, wireless connection (option)**
Display	5.7" colour, resistive touch screen	5.7" colour, resistive touch screen	5.7" colour, resistive touch screen	5.7" colour, resistive touch screen
Keypad	8 keys	8 keys	8 keys	8 keys
Protection class	IP 43	IP 43	IP 43	IP 43
Databases	19	19	19	19
Touch-free operation	2 programmable proximity sensors	2 programmable proximity sensors	2 programmable proximity sensors	2 programmable proximity sensors
USB-A	2	2	2	2
Ethernet	10 / 100 Mbit	10 / 100 Mbit	10 / 100 Mbit	10 / 100 Mbit
RS 232	2	2	2	2
Wi-Fi®	802.11 b/g/n	802.11 b/g/n	802.11 b/g/n	802.11 b/g/n
Power supply	13.5 ÷ 16 V DC	13.5 ÷ 16 V DC	13.5 ÷ 16 V DC	13.5 ÷ 16 V DC
Power consumption	10 W	10 W	10 W	10 W
Operating temperature	+10 ÷ +40 °C	+10 ÷ +40 °C	+10 ÷ +40 °C	+10 ÷ +40 °C
Operating temperature change rate	±0.3°C/1h (±1°C/8h)	±0.3°C/1h (±1°C/8h)	±0.3°C/1h (±1°C/8h)	±0.3°C/1h (±1°C/8h)
Atmospheric humidity***	40 ÷ 80%	40 ÷ 80%	40 ÷ 80%	40 ÷ 80%
Relative humidity change rate	±1%/h (±4%/8h)	±1%/h (±4%/8h)	±1%/h (±4%/8h)	±1%/h (±4%/8h)
Transport and storage temperature	-10 ÷ +50 °C	-10 ÷ +50 °C	-10 ÷ +50 °C	-10 ÷ +50 °C
Weighing pan dimensions	ø 16 mm	ø 16 mm	ø 60 mm (for filters), ø 16 mm	ø 16 mm
Weighing chamber dimensions	ø 90 × 90 mm	ø 90 × 90 mm	ø 90 × 90 mm	ø 90 × 90 mm
Weighing device dimensions	411 × 163 × 183 mm	411 × 163 × 183 mm	411 × 163 × 183 mm	411 × 163 × 183 mm
Net weight	9.1 kg	9.1 kg	9.1 kg	9.1 kg
Gross weight	16.6 kg	16.6 kg	16.6 kg	16.6 kg
Packaging dimensions	660 × 660 × 455 mm	660 × 660 × 455 mm	660 × 660 × 455 mm	660 × 660 × 455 mm
Anti-Draft Chamber Net Weight	7.6 kg	7.6 kg	7.6 kg	7.6 kg
Anti-Draft Chamber Gross weight	10.9 kg	10.9 kg	10.9 kg	10.9 kg
Anti-Draft Chamber Packaging dimensions	695 × 485 × 385 mm	695 × 485 × 385 mm	695 × 485 × 385 mm	695 × 485 × 385 mm

* Repeatability is expressed as a standard deviation from 10 weighing cycles. Values of parameters provided in Technical Specifications table have been determined under stable laboratory conditions. Due to ambient conditions impact or/and balance setup, the above parameters may vary for environments other than laboratory.

** optional solution on purchase order

*** non-condensing conditions

Wi-Fi® is a registered trademark of Wi-Fi® Alliance.

Technical Specifications

	MYA 5.4Y PLUS	MYA 6.4Y PLUS	MYA 11.4Y PLUS	MYA 11/52.4Y PLUS
Maximum capacity [Max]	5.1 g	6.1 g	11 g	11 g / 52 g
Minimum load	100 µg	100 µg	100 µg	100 µg
Readability [d]	1 µg	1 µg	1 µg	1 µg / 10 µg
Verification scale interval [e]	1 mg	1 mg	1 mg	1 mg
Tare range	-5.1 g	-6.1 g	-11 g	-52 g
Standard values				
Repeatability (5% Max)*	0.6 µg	0.6 µg	0.9 µg	1.5 µg
Repeatability (Max) *	1.6 µg	1.6 µg	2.5 µg	10 µg
Minimum weight (USP)	1.2 mg	1.2 mg	1.8 mg	3 mg
Minimum weight (U=1%, k=2)	0.12 mg	0.12 mg	0.18 mg	0.3 mg
Permissible values				
Repeatability (5% Max)*	1.2 µg	1.2 µg	1.6 µg	2.7 µg
Repeatability (Max)*	2.4 µg	2.4 µg	3.5 µg	15 µg
Linearity	±5 µg	±5 µg	±6 µg	±10 µg / ±30 µg
Eccentric load deviation	5 µg	5 µg	6 µg	6 µg / 10 µg
Stabilization time	max 8 s	max 8 s	max 10 s	max 10 s
Adjustment	internal	internal	internal	internal
Moveable range	—	—	—	Yes
Verification	Yes	Yes	Yes	Yes
OIML Class	I	I	I	I
Indicator fastening	35 cm cable, wireless connection (option)**	35 cm cable, wireless connection (option)**	35 cm cable, wireless connection (option)**	35 cm cable, wireless connection (option)**
Display	5.7" colour, resistive touch screen	5.7" colour, resistive touch screen	5.7" colour, resistive touch screen	5.7" colour, resistive touch screen
Keypad	8 keys	8 keys	8 keys	8 keys
Protection class	IP 43	IP 43	IP 43	IP 43
Databases	19	19	19	19
Touch-free operation	2 programmable proximity sensors	2 programmable proximity sensors	2 programmable proximity sensors	2 programmable proximity sensors
USB-A	2	2	2	2
Ethernet	10 / 100 Mbit	10 / 100 Mbit	10 / 100 Mbit	10 / 100 Mbit
RS 232	2	2	2	2
Wi-Fi®	802.11 b/g/n	802.11 b/g/n	802.11 b/g/n	802.11 b/g/n
Power supply	13.5 ÷ 16 V DC	13.5 ÷ 16 V DC	13.5 ÷ 16 V DC	13.5 ÷ 16 V DC
Power consumption	10 W	10 W	10 W	10 W
Operating temperature	+10 ÷ +40 °C	+10 ÷ +40 °C	+10 ÷ +40 °C	+10 ÷ +40 °C
Operating temperature change rate	±0.3°C/1h (±1°C/8h)	±0.3°C/1h (±1°C/8h)	±0.3°C/1h (±1°C/8h)	±0.3°C/1h (±1°C/8h)
Atmospheric humidity***	40 ÷ 80%	40 ÷ 80%	40 ÷ 80%	40 ÷ 80%
Relative humidity change rate	±1%/h (±4%/8h)	±1%/h (±4%/8h)	±1%/h (±4%/8h)	±1%/h (±4%/8h)
Transport and storage temperature	-10 ÷ +50 °C	-10 ÷ +50 °C	-10 ÷ +50 °C	-10 ÷ +50 °C
Weighing pan dimensions	ø 26 mm	ø 26 mm	ø 26 mm	ø 40 mm, ø 26 mm
Weighing chamber dimensions	ø 90 × 90 mm	ø 90 × 90 mm	ø 90 × 90 mm	ø 90 × 90 mm
Weighing device dimensions	411 × 163 × 183 mm	411 × 163 × 183 mm	411 × 163 × 183 mm	411 × 163 × 183 mm
Net weight	9.1 kg	9.1 kg	9.1 kg	9.1 kg
Gross weight	16.6 kg	16.6 kg	16.6 kg	16.6 kg
Packaging dimensions	660 × 660 × 455 mm	660 × 660 × 455 mm	660 × 660 × 455 mm	660 × 660 × 455 mm
Anti-Draft Chamber Net Weight	7.6 kg	7.6 kg	7.6 kg	7.6 kg
Anti-Draft Chamber Gross weight	10.9 kg	10.9 kg	10.9 kg	10.9 kg
Anti-Draft Chamber Packaging dimensions	695 × 485 × 385 mm	695 × 485 × 385 mm	695 × 485 × 385 mm	695 × 485 × 385 mm

* Repeatability is expressed as a standard deviation from 10 weighing cycles. Values of parameters provided in Technical Specifications table have been determined under stable laboratory conditions. Due to ambient conditions impact or/and balance setup, the above parameters may vary for environments other than laboratory.

** optional solution on purchase order

*** non-condensing conditions

Wi-Fi® is a registered trademark of Wi-Fi® Alliance.

Technical Specifications

	MYA 21.4Y PLUS	MYA 31.4Y PLUS	MYA 21/52.4Y PLUS
Maximum capacity [Max]	21 g	31 g	21 g / 52 g
Minimum load	100 µg	100 µg	100 µg
Readability [d]	1 µg	1 µg	1 µg / 10 µg
Verification scale interval [e]	1 mg	1 mg	1 mg
Tare range	-21 g	-31 g	-52 g
Standard values			
Repeatability (5% Max)*	1 µg	1.2 µg	1.5 µg
Repeatability (Max)*	3 µg	4.5 µg	10 µg
Minimum weight (USP)	2 mg	2.4 mg	3 mg
Minimum weight (U=1%, k=2)	0.2 mg	0.24 mg	0.3 mg
Permissible values			
Repeatability (5% Max)*	1.6 µg	2.5 µg	2.7 µg
Repeatability (Max)*	4 µg	6.5 µg	15 µg
Linearity	±7 µg	±8 µg	±10 µg / ±30 µg
Eccentric load deviation	7 µg	8 µg	6 µg / 10 µg
Stabilization time	max 10 s	max 10 s	max 10 s
Adjustment	internal	internal	internal
Moveable range	—	—	—
Verification	Yes	Yes	Yes
OIML Class	I	I	I
Indicator fastening	35 cm cable, wireless connection (option)**	35 cm cable, wireless connection (option)**	35 cm cable, wireless connection (option)**
Display	5.7" colour, resistive touch screen	5.7" colour, resistive touch screen	5.7" colour, resistive touch screen
Keypad	8 keys	8 keys	8 keys
Protection class	IP 43	IP 43	IP 43
Databases	19	19	19
Touch-free operation	2 programmable proximity sensors	2 programmable proximity sensors	2 programmable proximity sensors
USB-A	2	2	2
Ethernet	10 / 100 Mbit	10 / 100 Mbit	10 / 100 Mbit
RS 232	2	2	2
Wi-Fi®	802.11 b/g/n	802.11 b/g/n	802.11 b/g/n
Power supply	13.5 ÷ 16 V DC	13.5 ÷ 16 V DC	13.5 ÷ 16 V DC
Power consumption	10 W	10 W	10 W
Operating temperature	+10 ÷ +40 °C	+10 ÷ +40 °C	+10 ÷ +40 °C
Operating temperature change rate	±0.3°C/1h (±1°C/8h)	±0.3°C/1h (±1°C/8h)	±0.3°C/1h (±1°C/8h)
Atmospheric humidity***	40 ÷ 80%	40 ÷ 80%	40 ÷ 80%
Relative humidity change rate	±1%/h (±4%/8h)	±1%/h (±4%/8h)	±1%/h (±4%/8h)
Transport and storage temperature	-10 ÷ +50 °C	-10 ÷ +50 °C	-10 ÷ +50 °C
Weighing pan dimensions	ø 26 mm	ø 26 mm	ø 40 mm, ø 26 mm
Weighing chamber dimensions	ø 90 × 90 mm	ø 90 × 90 mm	ø 90 × 90 mm
Weighing device dimensions	411 × 163 × 183 mm	411 × 163 × 183 mm	411 × 163 × 183 mm
Net weight	9.1 kg	9.1 kg	9.1 kg
Gross weight	16.6 kg	16.6 kg	16.6 kg
Packaging dimensions	660 × 660 × 455 mm	660 × 660 × 455 mm	660 × 660 × 455 mm
Anti-Draft Chamber Net Weight	7.6 kg	7.6 kg	7.6 kg
Anti-Draft Chamber Gross weight	10.9 kg	10.9 kg	10.9 kg
Anti-Draft Chamber Packaging dimensions	695 × 485 × 385 mm	695 × 485 × 385 mm	695 × 485 × 385 mm


* Repeatability is expressed as a standard deviation from 10 weighing cycles. Values of parameters provided in Technical Specifications table have been determined under stable laboratory conditions. Due to ambient conditions impact on balance setup, the above parameters may vary for environments other than laboratory.

** optional solution on purchase order

*** non-condensing conditions

Wi-Fi® is a registered trademark of Wi-Fi® Alliance.

Dimensions

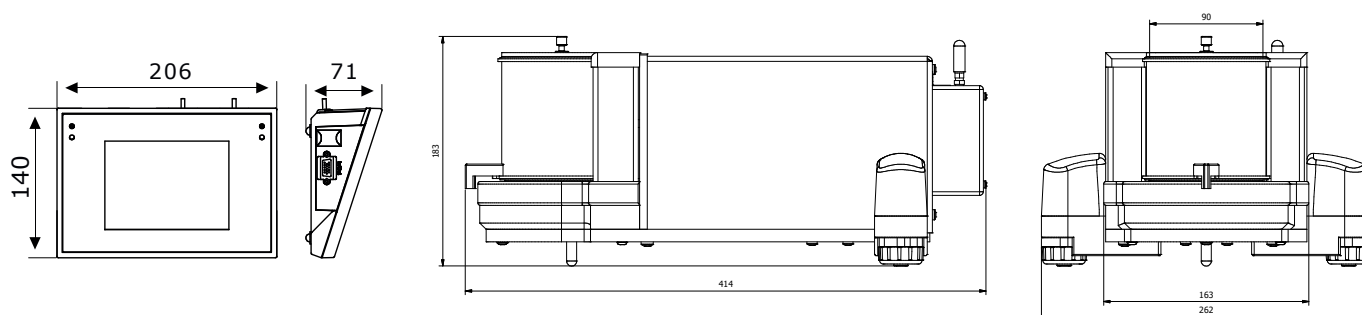

Accessories

Weighing Tables

- granite antivibration table
- antivibration tables for laboratory balances

Professional weighing

- Adapter for calibration of MY11 series pipettes

Ambient Conditions

- DJ-05 anti-static ionizer
- THB-S or THB-P sensor

Peripheral Devices

- Epson dot matrix printer

- barcode scanners

- WD-5/3Y LCD display

Cables, Converters

- P0108: RS 232 cable (balance-computer)
- P0167: RS 232 cable (balance-computer)
- P0151: RS 232 cable (balance - Epson printer)

Electrical Accessories

- ZR-02 power supply with battery

Dedicated Software

R-LAB

- collecting measurements
- carrying out statistical analysis of measurements
- customized graphs and reports

E2R Weighing Records

- complete, automated databases synchronization
- fully supported processes of labelling and parts counting
- record of weighings, weighings archiving
- basic and advanced (with graphs) reports

Label Editor R02

- designing label templates
- sending graphics and fonts to label printers
- printing label templates using connected printers

Pipettes

- determining measurement errors of pipettes volume
- accordance with ISO 8655
- calibration of single-channel and multi-channel pipettes
- calibration of fixed-volume and variable-volume pipettes

Audit Trail Reader

- support of Audit Trail function available for 3Y, 4Y, HY10, WLY, WPY series weighing instruments
- record of operator's activity from the moment of logging in

Parameters Editor

- remote change of parameters
- remote on-line preview of the display
- displaying current mass indication
- software update
- file loading, editing and saving parameters to a file
- import and export of parameters
- interfaces: RS232, Ethernet and Wireless Connection.
- quick and easy edition of balance parameters using computer.

RAD KEY

- Establishing cooperation between a weighing instrument and a computer

R. Barcode

- The basic function software is presentation of the data sent by barcode scanners connected to PC via USB or RS232

Radwag Development Studio

- presentation of functions (and subfunctions) of communication protocol (Common Communication Protocol)
- possibility of connection with weighing equipment on which each function is carried out,
- library with mass control, contained within the development environment
- complete documentation of the communication protocol
- set of user manuals for different solutions addressed for programmers employed in companies using RADWAG-manufactured weighing equipment

LabView Driver

- operation of RADWAG balances in LabView environment

RADWAG Connect

- establishing communication with all balances, scales and weighing modules using Common Communication Protocol
- communication via local network,
- support of basic functions
- auto searching for devices
- connecting with few devices simultaneously, swapping between them
- clear list of connected platforms
- record of measurements in the program,
- export of carried out measurements to CSV file,
- work performed using freely selected device with Windows 10 operating system

RADWAG Remote Desktop

- remote operation via computer, mobile phone or tablet
- sending text messages
- version for Windows 10 and Android systems