

DRY CABINETS

Vacuum Desiccators



4706/0063, Polycarbonate Vacuum Desiccator

- Main unit and lid are completely transparent. You can instantly know the status of items that are stored and when replacing silica gel.
- Dry seal type uses an O-ring. Eliminated the trouble of applying grease.
- Maximum strength in the transparent resins. Difficult to crack even if dropped, so you can use it with confidence.
- The ceramic plate is heat resistant, corrosion resistant, has an anti-sticking property and won't scratch.
- As a feature of PC, it has strong tensile strength, is durable to impacts, & is heat-resistant.


0063 - With Spigot, Effective inside

Diameter: 230mm

4706 - With Spigot, Vacuum

Gauge, Effective inside Diameter:

230mm

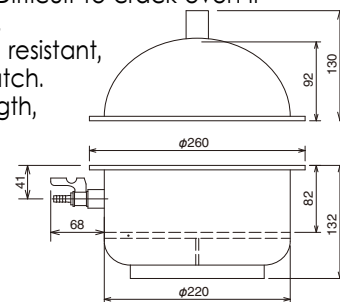

Model	4706	0063
Body inside diameters	210mm	250mm
Body inside height	90mm H	100mm H
A dedicated ceramic plate	O	O
Spigot nozzle diameter	6 to 8 mm	6 to 8 mm
Spigot material	Polycarbonate	Polycarbonate
Vacuum gauge	0 to -0.1 MPa (minimum scale units: 0.005 MPa)	-